

Faculty of Engineering and Technology


Applied Chemistry Laboratory
Subject: Engineering Chemistry

Observation:

PART 1 - WELL WATER

Screen shot-

Water Analysis : Chemical Content

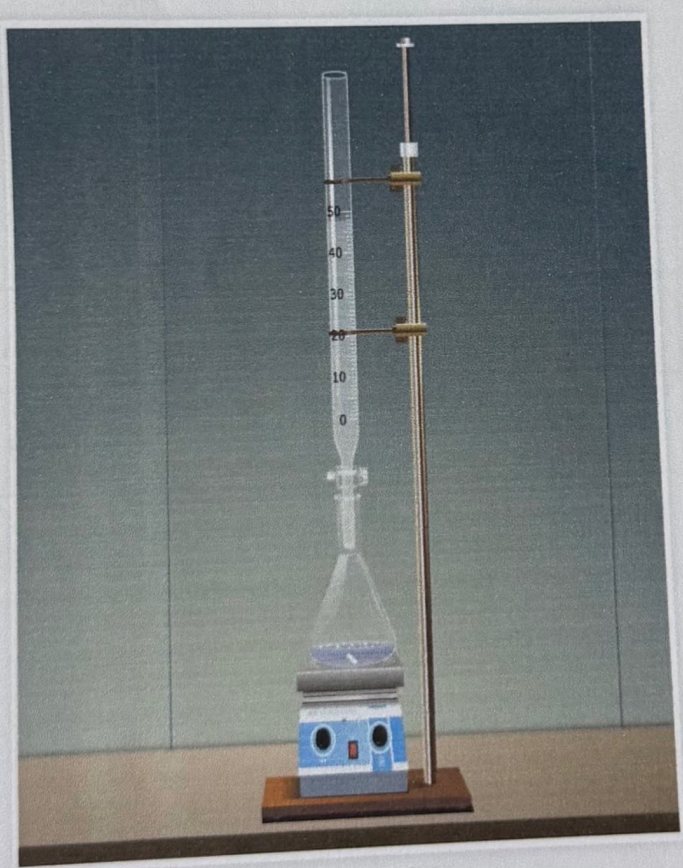

Copyright © Amrita University 2005 - 2015

VARIABLES

Select test:
Hardness

Titrate :
Well water

Volume of titrate : 10

Speed of titrant : 0.2

Molarity of titrant : 0.01

START

RESET

RESULT

Titrant : EDTA
Titrant used : 2.5 ml
Indicator : Eriochrome Black T

Burette

Conical flask

Indicator



End point

: 0.01 M EDTA solution
: 10 mL of sample + Indicator
: EBT
: Wine red to Blue

Faculty of Engineering and Technology

Applied Chemistry Laboratory
Subject: Engineering Chemistry

Reaction:

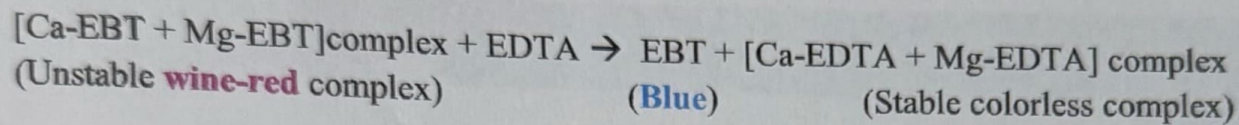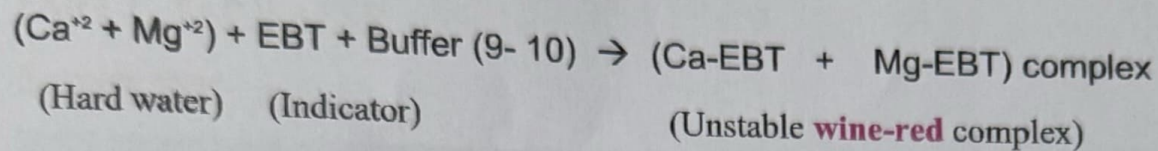

Pilot Reading : 2(mL) to 3(mL)

Reading	I (mL)	II (mL)	Constant (mL)
Initial	0.00	0.00	2.5
Final	2.5	2.5	
Difference	2.5	2.5	


Faculty of Engineering and Technology

Applied Chemistry Laboratory
Subject: Engineering Chemistry

PART 2 - TAP WATER

Screen shot-

Water Analysis : Chemical Content

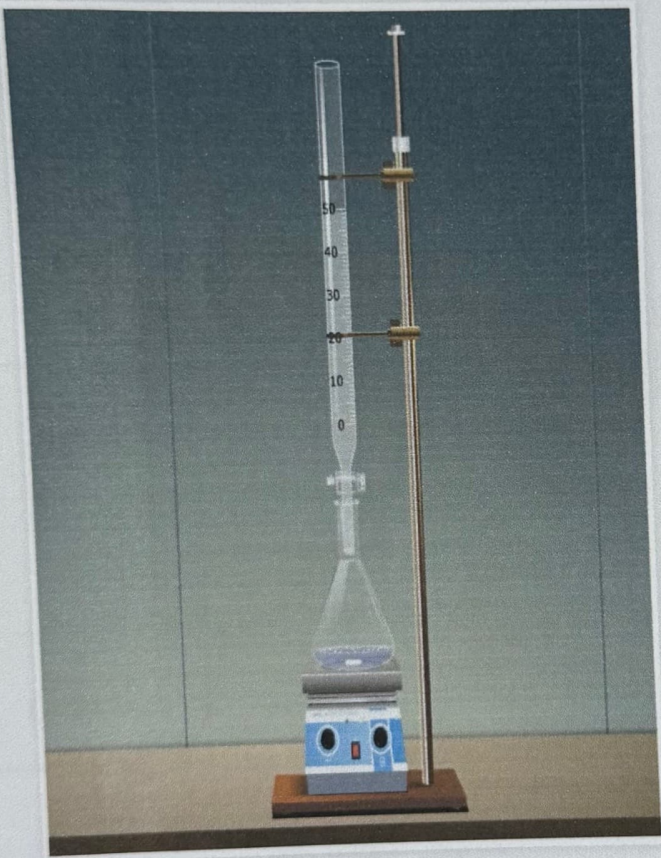

VARIABLES

Select test:
Hardness

Titrate:
Tap water

Volume of titrate : 12

Speed of titrant : 0.4

Molarity of titrant : 0.02

START

RESET

RESULT

Titrant : EDTA

Titrant used : 0.9 ml

Indicator : Eriochrome Black T

Copyright © Azimkha University 2005 - 2015

Burette

Conical flask

Indicator

End point

: **0.02 M** EDTA solution

: **12 mL** of sample + Indicator

: EBT

: Wine red to Blue

Faculty of Engineering and Technology**Applied Chemistry Laboratory**
Subject: Engineering Chemistry**Reaction:**

$\text{Ca}^{+2} + \text{Mg}^{+2} + \text{EBT} + \text{Buffer (pH 9-10)} \rightarrow (\text{Ca-EBT} + \text{Mg-EBT}) \text{ complex}$
(Ca (Hard water) (Indicator) (Unstable **wine-red** complex)

$[\text{Ca-EBT} + \text{Mg-EBT}] \text{ complex} + \text{EDTA} \rightarrow \text{EBT} + [\text{Ca-EDTA} + \text{Mg-EDTA}] \text{ complex}$
(Unstable **wine-red** complex) (Blue) (Stable colourless complex)

Pilot Reading : 0.5 (mL) to 1 (mL)

Reading	I (mL)	II (mL)	Constant (mL)
Initial	0.00	0.00	<div>1.9</div>
Final	1.9	1.9	
Difference	1.9	1.9	


Faculty of Engineering and Technology

Applied Chemistry Laboratory
Subject: Engineering Chemistry

PART 3 - SEA WATER

Screen shot-

Water Analysis : Chemical Content

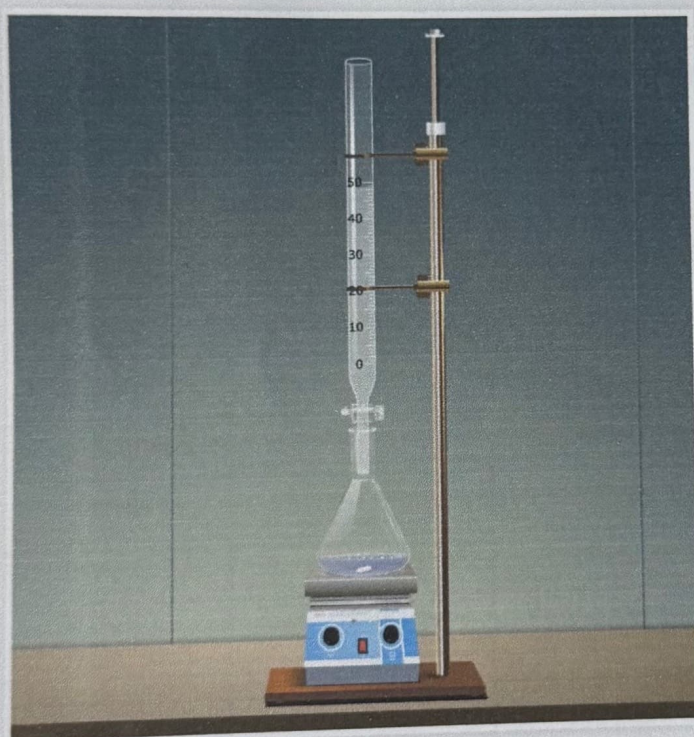

VARIABLES

Select test:

Hardness

Titrate:

Sea water

Volume of titrate : 13

Speed of titrant : 0.5

Molarity of titrant : 0.05

START

RESET

RESULT

Titrant : EDTA

Titrant used : 1.6 ml

Indicator : Eriochrome Black T



Burette
Conical flask
Indicator
End point

: 0.05M EDTA solution
: 13mL of sample + Indicator
: EBT
: Wine red to Blue

Faculty of Engineering and Technology

Applied Chemistry Laboratory
Subject: Engineering Chemistry

Reaction:

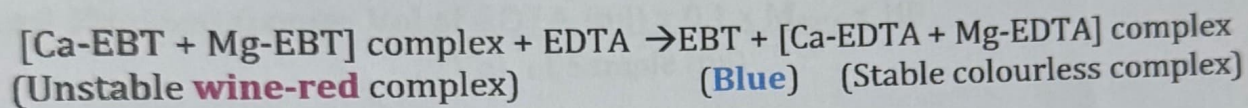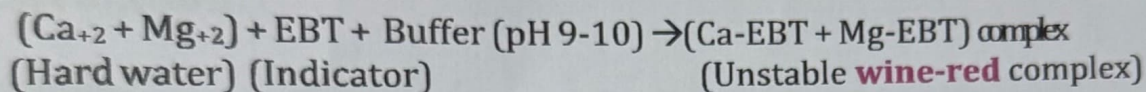

Pilot Reading : 1 (mL) to 2 (mL)

Reading	I (mL)	II (mL)	Constant (mL)
Initial	0.00	0.00	1.6
Final	1.6	1.6	1.6
Difference	1.6	1.6	1.6

Faculty of Engineering and Technology

Applied Chemistry Laboratory
Subject: Engineering Chemistry

Calculation:

Formula:

$$\text{Total Hardness (ppm)} = \frac{\text{Vol. of EDTA (ml)} \times 0.1 \times M_{\text{EDTA}} \times 10^6}{\text{Vol. of Sample (ml)}}$$

Part-1: Well water

$$\begin{aligned} \text{Total Hardness (ppm)} &= \frac{2.5 \times 0.1 \times 0.01 \times 10^6}{10} \\ &= \underline{250} \text{ ppm} \end{aligned}$$

Part-2: Tap water

$$\begin{aligned} \text{Total Hardness (ppm)} &= \frac{0.9 \times 0.02 \times 0.1 \times 10^6}{12} \\ &= \underline{150} \text{ ppm} \end{aligned}$$

Part-3: Sea water

$$\begin{aligned} \text{Total Hardness (ppm)} &= \frac{1.6 \times 0.1 \times 0.05 \times 10^6}{13} \\ &= \underline{615.38} \text{ ppm} \end{aligned}$$

Result:

The hardness of Sea water is the highest = 615.38 ppm.

Faculty of Engineering and Technology**Applied Chemistry Laboratory**
Subject: Engineering Chemistry**Assignments:**

1. Convert the total hardness of water samples in ppm and meq/L from the following:
A. 20.23_{CaCl_2} ; B. 31.8_{Fe}
2. A sample of water has hardness 208ppm CaCO_3 eq. Find the hardness in terms of mg/L, CaCl_2 , meq/L.
3. How many grams of FeSO_4 dissolved per litre gives 210.5ppm of hardness?
4. Three samples were analysed for their salt content:
 - (i) Sample A contains 168mg of magnesium carbonate per litre
 - (ii) Sample B contains 820mg of calcium nitrate and 2mg of silica per litre
 - (iii) Sample C contains 20g potassium nitrate and 2g calcium carbonate per 500mlDetermine the hardness in ppm and grains per gallon.
5. Classify the following into temporary, permanent and non-hardness causing impurities: $\text{Ca}(\text{HCO}_3)_2$, MgSO_4 , CaCl_2 , CO_2 , HCl , $\text{Mg}(\text{HCO}_3)_2$, CaSO_4 , NaCl .
How many grams of CaCl_2 dissolved per litre gives 150ppm of hardness?
6. Classify the following into carbonate and non-carbonate impurities and calculate all types of hardness-
 $\text{Mg}(\text{HCO}_3)_2 = 7.1\text{mg/L}$, $\text{Ca}(\text{HCO}_3)_2 = 8.1\text{mg/L}$, $\text{MgCO}_3 = 4.2\text{mg/L}$, $\text{CaCO}_3 = 10\text{mg/L}$,
 $\text{MgSO}_4 = 24\text{mg/L}$.
7. 0.28g CaCO_3 was dissolved in HCl and made upto 1L with distilled water. 100ml of this solution required 28ml EDTA solution. 100ml of hard water sample required 33ml of EDTA solution. After boiling, cooling and filtering, 100ml of the solution required 10ml of EDTA. Calculate hardness.